


Initial Fitting - for Medical Professionals:



1. Open up the knee brace. Set the desired length of the stays with help of the telescopic mechanism. Ensure that the stays are of equal length on both sides. If necessary, the leg cuffs can be moved.

Instruction video:

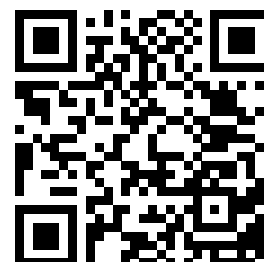


2. Bend the aluminum stays to match the contour of the leg. Place the joint on a stable surface and bend the stays carefully but firmly. **Note!** Move the aluminum cuffs before bending the stays. Then slide them back to the desired position.



3. Align the orthosis hinge with the knee joint. Then fasten the strap buckles in the order shown (A-D) and tighten the straps. The strap may be cut to desired length.

Make sure the patient is not experiencing any discomfort.

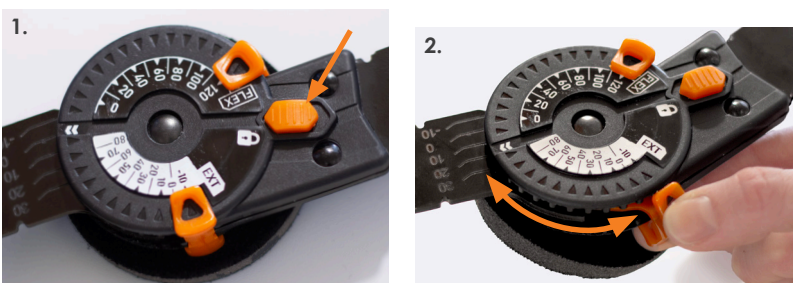

ROM-settings:

SE 1. The joint can be locked with the quick-release button from -10° to 30° at the time of fitting. When the button is unlocked, the patient can move the joint within the set range of motion. Repeat the adjustments on both sides.

2. Set the desired range of motion in the joint; pull the adjustment button and slide to the desired position, release the button. Repeat the setting on both sides.