

SOLUTIONS

bite registration

HARD-HITTING SOLUTIONS !



Clinical

Bite Registration Solutions

Zhermack®

CLINICAL



bite registration solutions

Imperceptible consistency with extremely rigid setting.

ZHERMACK presents a range of *bite registration solutions* with exclusive features.

occlufast rock

colorbite rock

colorbite D

bite registration



IMPERCEPTIBLE CONSISTENCY

› When patients close their mouth they do not have *any perception of the material's consistency*.

RIGIDITY

› The optimal final rigid set *prevents compressions*, incorrect re-positioning and allows *easy trimming*.

SPEED

› Minimal time in the mouth ensures *greater patient comfort* and *reduces possible errors* due to minor jaw movements.

THIXOTROPIC

› **ZHERMACK** bite registration solutions allow *easy control* and positioning.

colorbite D – colorbite rock

The innovation that changes the bite registration.

bite registration

COLORBITE D

Extra rigid, thermochromic A-silicone with fast set.
COLORBITE D is **green** at room temperature and becomes **yellow** at 35°C in the mouth to indicate that setting is complete.



Application of **COLORBITE D** (green)



Setting phase (colour change from green to yellow)



Complete colour change (yellow)

FEATURES

- **Imperceptible** consistency
- **Thixotropic**, does not drip
- Extremely **rigid after setting** (45 Shore D)
- **Fast set**
- Very **precise**
- **Easy to finish**
- **Lime flavour**

ADVANTAGES

- The use of thermochromic pigments allows you **see** if the material has set, therefore minimising the margin of error
- **Imperceptible to the patient** to give a natural bite
- After setting it can be **finished** using a knife or be milled
- **Easy to position** both in occlusal and buccal techniques

TECHNICAL DATA

Colour change	green -----> yellow
Time in the mouth	1'
Reproduction of details	2 µm
Dimensional variation (after 24 hours)	< -0,05 %
Final hardness	45 Shore D

COLORBITE ROCK

Rigid, thermochromic A-silicone with extra fast set.
COLORBITE ROCK is **purple** at room temperature and becomes **pink** at 35°C in the mouth to indicate that setting is complete.



Application of **COLORBITE ROCK** (purple)



Setting phase (colour change from purple to pink)



Complete colour change (pink)

FEATURES

- Extra **quick**
- **Imperceptible** consistency
- **Thixotropic**, does not drip
- **Rigid after setting** (95 Shore A)
- Very **precise**
- **Easy to finish**
- **Raspberry flavour**

ADVANTAGES

- The use of thermochromic pigments allows you **see** if the material has set, therefore minimising the margin of error
- **Imperceptible to the patient** to give a natural bite
- After setting it can be **finished** using a knife or be milled
- **Easy to position** both in occlusal and buccal techniques

TECHNICAL DATA


Colour change	purple -----> pink
Time in the mouth	45"
Reproduction of details	2 µm
Dimensional variation (after 24 hours)	< -0,05 %
Final hardness	95 Shore A

occlufast rock

The Number 1 bite registration solution in the world.

FEATURES

- **Imperceptible consistency:** when patients close their mouth they do not have *any perception of the material's consistency*
- **95 Shore A rigidity:** its rigidity prevents compression and incorrect re-positioning. Allows easy trimming or milling
- **Fast 60 secs:** minimal time in the mouth ensures *greater patient comfort* and *reduces possible errors* due to minor jaw movements
- **Thixotropic:** **OCCLUFAST ROCK** allows easy control and positioning

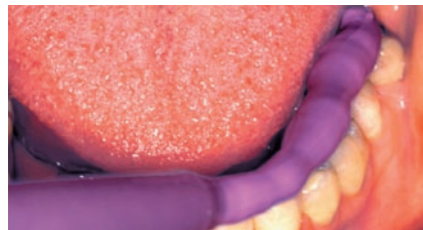

Application of **OCCLUFAST ROCK** (first phase)



Bite registration (second phase)



Buccal registration

ADVANTAGES

- **Easy to use** and **no air bubbles**
- **Control** of positioning, does not drip
- Minimum thickness, virtually equal to zero
- **Easy to trim or mill** to eliminate excess material

TECHNICAL DATA

Time in the mouth	1'
Reproduction of details	2 µm
Dimensional variation (after 24 hours)	< -0,05 %
Final hardness	95 Shore A



bite registration solutions
Rigidity in the pure state!

bite registration



Packaging

C200726 OCCLUFAS ^T ROCK	2x50 ml cartridges + 12 mixing tips (green)
C200750 COLORBITE D	2x50 ml cartridges + 12 mixing tips (green)
C200760 COLORBITE ROCK	2x50 ml cartridges + 12 mixing tips (green)

Accessories

C202100 DISPENSER D2	D2 Dispenser - 1:1
U100210 PLUNGER	Plunger for Dispenser D2 - 1:1
C202080 MIXING TIPS STANDARD	Mixing tips, green (48 pcs)
C202095 BITE TIPS	Intra-oral bite tips (48 pcs)

Branch
Zhermack GmbH Deutschland
 Öhlmühle 10 - D-49448 Marl
 Tel. +49 - 0 54 43 / 20 33 - 0
 Fax +49 - 0 54 43 / 20 33 - 11
 info@zhermack.de
 www.zhermack.com

Branch
Zhermapol
 ul. Augustówka 14
 02 - 981 Warszawa - Polska
 Tel. +48 - 22 858 82 72 - 858 73 41
 Fax +48 - 22 642 07 14
 zhermapol@zhermapol.pl
 www.zhermack.com

Branch
Zhermack Inc.
 P.O. Box 4195
 River Edge, New Jersey 07661-4195
 Tel. +1 (877) 819-6206
 +1 (732) 389 8540
 Fax +1 (732) 389 8543
 info@zhermackusa.com
 www.zhermackusa.com

Representative Office
Zhermack UK and Ireland
 21 Balmoral Road
 Melton Mowbray - Leics.
 LE13 1JN
 Tel. +44 - (0)7870 690811
 uk@zhermack.com
 www.zhermack.com

Representative Office
Zhermack Russia
 5/2 Marshala Vasilevskogo St.
 123098 Moscow - Russia
 Tel. +7 - 0916 9301191
 info.russia@zhermack.com
 www.zhermack.com

Zhermack SpA
 Via Bovazecchino, 100 - 45021 Badia Polesine (RO) ITALY
 Tel. +39 - 0425 597611 Fax +39 - 0425 53596 info@zhermack.com www.zhermack.com

