

Model 3200

# Pool lift

Direct from wheelchair to pool



- Direct transfers to and from wheelchair to water.
- Easy to move and store the lift.
- Use one lift at multiple locations: multiple pools or several areas at same pool.
- Variety of accessories make transfers easy regardless of disability.
- Designed for the swimming pool environment.

Lift & Care systems for people with reduced mobility

[www.handimove.com](http://www.handimove.com)

## Pool lift

### DIRECT TRANSFERS TO AND FROM WHEELCHAIR TO WATER.

The Pool lift brings the user to the water in one smooth motion. This efficient lift moves the user to the pool in one transfer.



The Pool lift has a pivoting range of 360°.

### EASY TO MOVE AND STORE THE LIFT.

Operation is simple, with the push of a button, the spindle lowers into the floor socket at the edge of the pool.

When finished, lift the spindle out of the socket, and wheel the lift to storage.

Once removed, seal the socket opening with the plug.

### USE ONE LIFT AT MULTIPLE LOCATIONS: MULTIPLE POOLS OR SEVERAL AREAS AT SAME POOL.

The installation of additional floor sockets, and the fact that it is easily movable, allows the lift to be used in different locations

### VARIETY OF ACCESSORIES MAKE TRANSFERS EASY REGARDLESS OF DISABILITY.

Customize the transfer with the right accessory.


Based on the need choose a pool lift chair, Handi-Move Body Support, bathing sling or standard sling.

### DESIGNED FOR THE SWIMMING POOL ENVIRONMENT.

Made of stainless steel and powered by batteries, the mobile Pool lift is designed for safe, long-term use in the pool environment.



DISCRETE SOCKET,  
WITH PLUG TO CAP SOCKET WHEN NOT IN USE



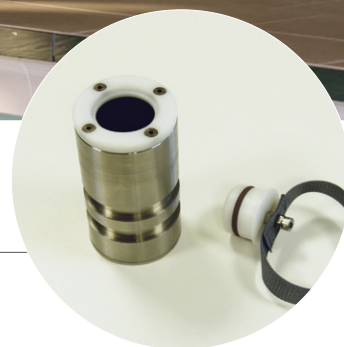





# Pool lift

## DIMENSIONS AND SPECIFICATIONS

| Dimensions (mm) |      | Specifications      |                     |
|-----------------|------|---------------------|---------------------|
| A               | 1840 | Safe working load   | 135 kg              |
| B               | 1305 | Batteries           | 2 x 12 Vdc; 8 Ah    |
| C               | 540  | Charger             | separate            |
| D               | 350  | Diameter of casters | 100 mm              |
| E               | 1700 | Construction        | stainless steel 316 |
| F               | 550  | Hand control        | electrical (IP67)   |
| G               | 500  |                     |                     |
| H               | 310  |                     |                     |
| I               | 370  |                     |                     |

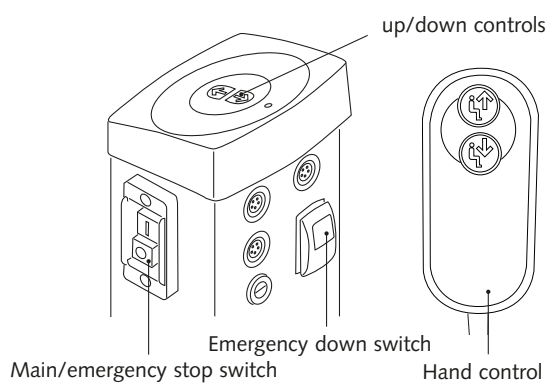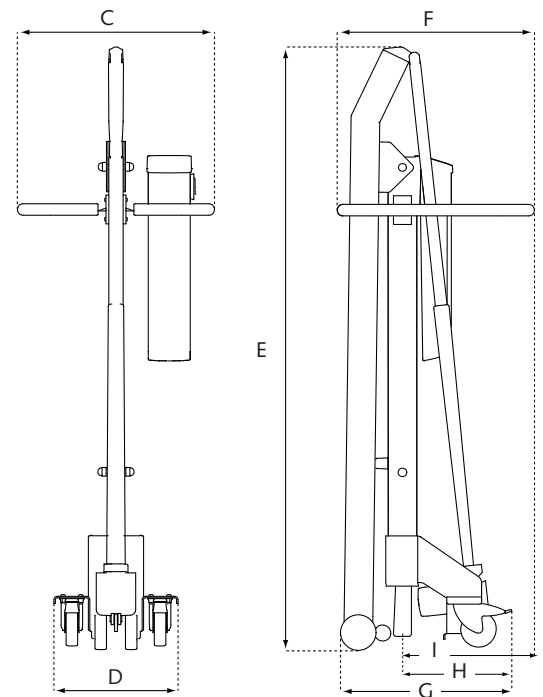
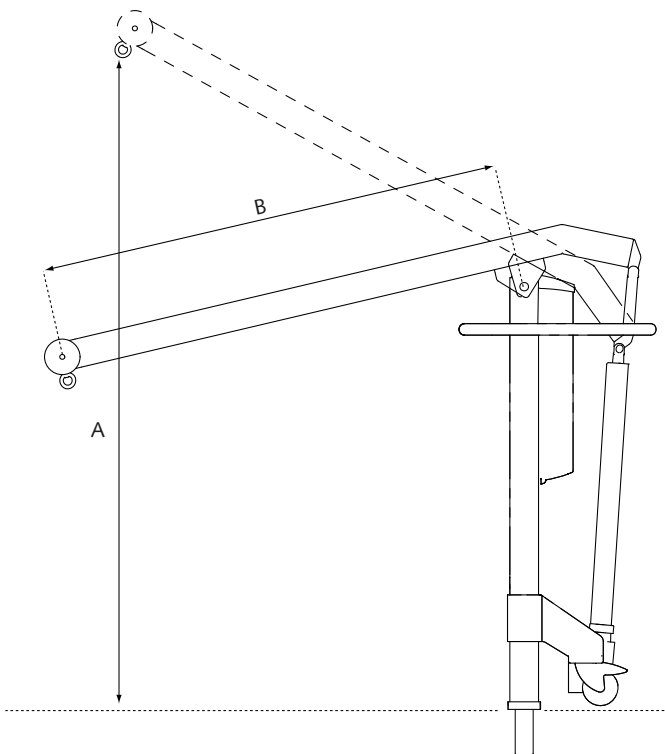


## ● BODY SUPPORT SYSTEMS

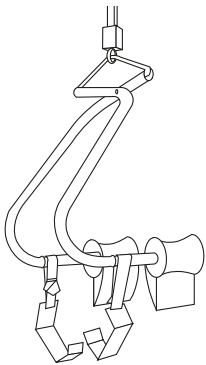



### Body Support

XS - S - M • art. 1035  
Electrolytically coated (anticorrosive)  
+ Slipproof Thigh Supports • art. 1074

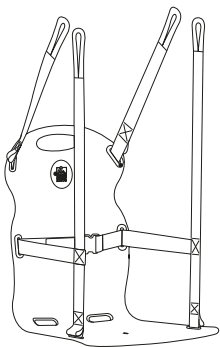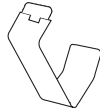

### Pool Lift Chair (PVC)



Uni-size • art. 3620  
+ Spreaderbar • art. 2110 / 2125

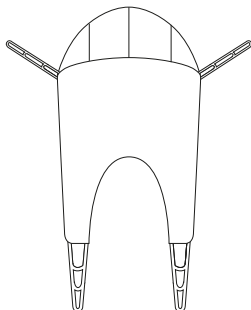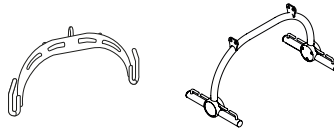

### Standard Sling


with moon shaped head support and padded thigh flaps

XS - S - M - L - XL • art. 3431  
+ Spreaderbar • art. 2110 / 2125

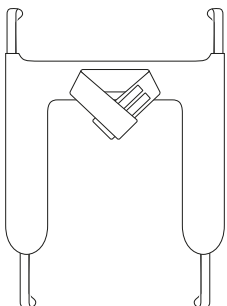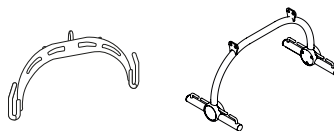

### Bathing Sling (PVC)

XS - S - M - L - XL • art. 3464  
+ Spreaderbar • art. 2110 / 2125

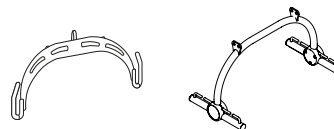

For free demonstration or further information

