

A woman with blonde hair, wearing a brown sweater and beige trousers, is walking on a city street. She is wearing a black Stratus stabilizing airwalker on her right foot. The background shows a modern city street with buildings and other pedestrians.

Stratus


Stabilizing airwalker

- ✓ Lightweight
- ✓ Slim design
- ✓ Low profile for a level posture
- ✓ Two built-in pumps
- ✓ Dual independent air chambers
- ✓ Even compression

NordiCare[®] 

Smooth walking with ankle immobilization

The Stratus Air Walker combines outstanding stability with all-day comfort to support recovery following foot and ankle injuries. Its lightweight, low-profile design promotes a more level and balanced walking pattern while providing reliable immobilization and protection throughout the treatment period.

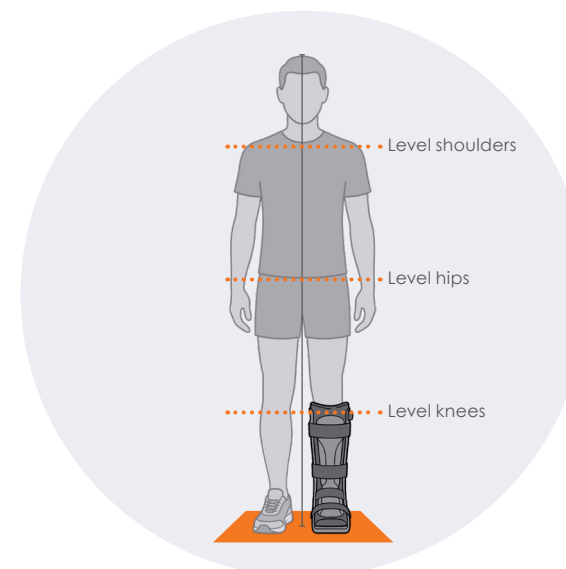

Low build height for level walking

The heel height align the walker with standard footwear, helping maintain a level and balanced gait and reducing compensatory tilting. This can help reduce strain on the back, hips and knees that may result from differences in leg height. The low profile reduces the need for compensation in the opposite shoe.

Individually adjustable fit

The fully integrated dual air pump system allows the fit and compression to be adjusted quickly and easily, without a separate pump. Two independent air chambers, one on each side, enable customized support tailored to the patient's needs, helping to maximize both comfort and stability. Four adjustable straps ensure a secure fit for a wide range of foot and leg shapes while the ankle strap plays a key role in maintaining proper foot positioning.



Effective compression

Evenly distributed pneumatic compression comfortably surrounds and immobilizes the foot and ankle, helping to reduce swelling while enhancing support and patient comfort throughout recovery.

Comfort and function in everyday life


The anatomically contoured shell and soft padding provide exceptional comfort during extended wear. Ventilated materials improve airflow (1) and the removable toe cover provides extra protection and warmth when needed (2). The liner is easy to remove and reapply. Included extra padding helps optimize fit and comfort where needed (3). The non-slip outsole enhances stability, traction and safety during everyday activities (4).



Additional Accessories

A range of accessories is available separately to further enhance the fit, comfort and functionality.

	Accessories	Size	Art.no
1	Socks, 2 pcs	Universal	802153000
2	Rain cover	S/M, L/XL	802200040 / 60
3	Hygiene protection, 10 pcs	Universal	802190000
4	Extra liner	XS-XL	802254030-70
5	Walker Wedge Kit 10°	1, 2	802240040 / 60

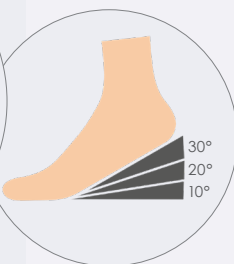

Adjustable Plantar Flexion with Walker Wedge Kit 10°

Easily adjust the Stratus Walker during Achilles tendon rehabilitation. The precisely angled, platform-free 10° wedges provide controlled positioning of the foot and ankle, helping to achieve the intended plantar flexion during recovery. Three wedges can be stacked to achieve 10°, 20° or 30° of plantar flexion.



Stratus Indications and art.no

- Stable foot and ankle fractures
- Severe ankle sprains and soft tissue injuries
- Postoperative recovery and immobilization
- Achilles tendon ruptures (with heel wedges)

Art.no	Size	Shoe size
802250030	XS	<35
802250040	S	36-39
802250050	M	40-44
802250060	L	44-46
802250070	XL	>46