

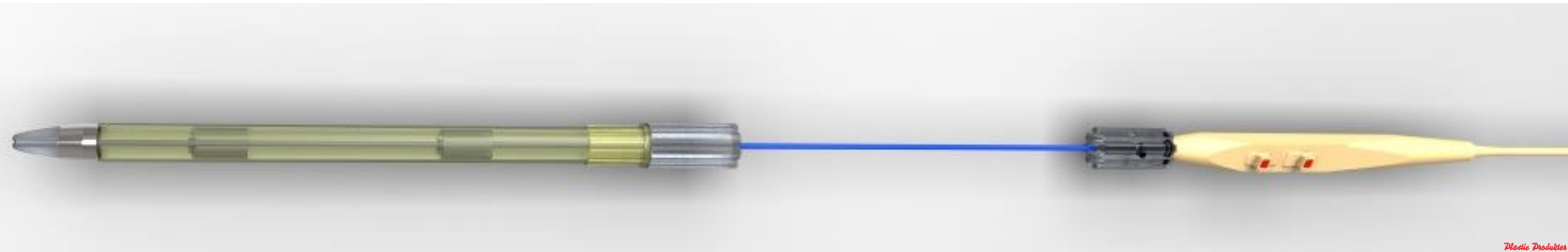
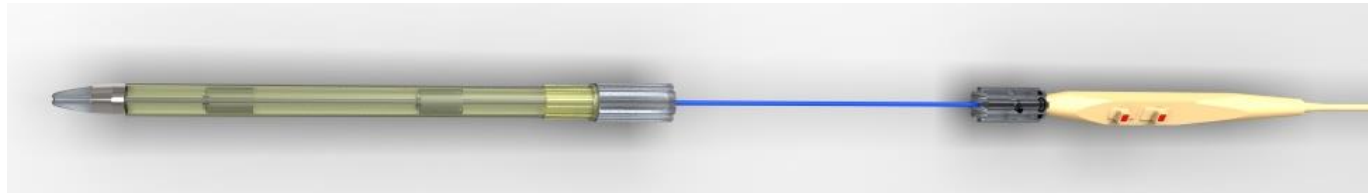
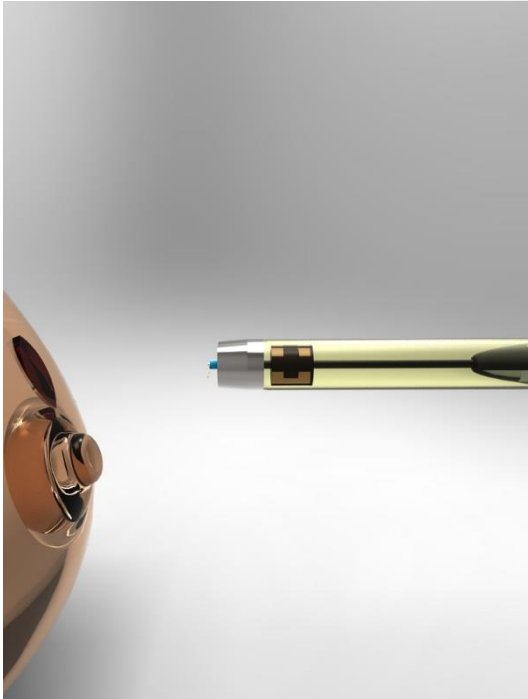
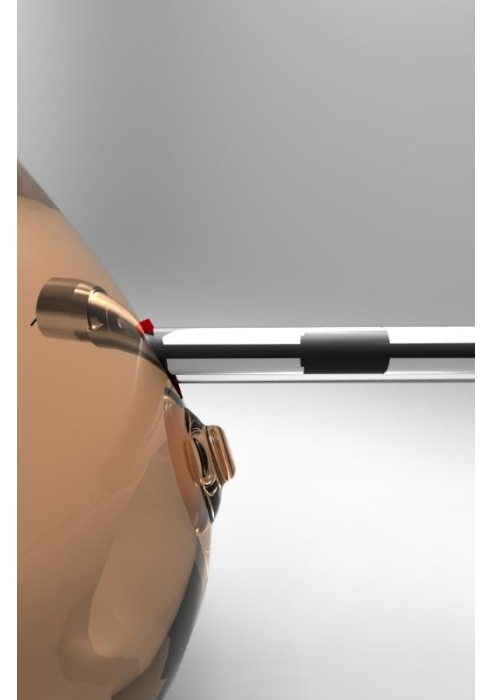


# NOVOAIM BIOPSY BY DIATERMI




Novoaim excitu kit with diathermy handle

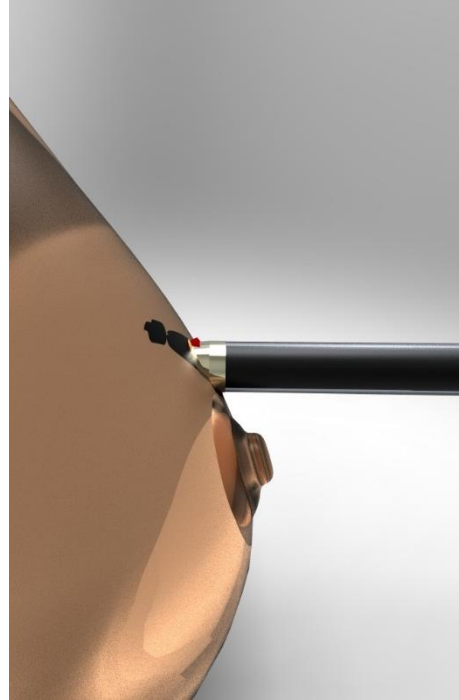


*Plastic Produkter*

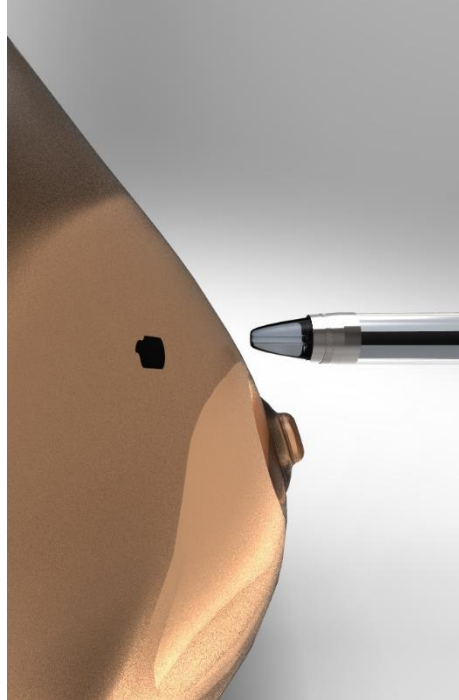


Disposable excitu biopsy kit, with diathermy.  
System and device development for NOVOAIM.  
The system is patented.

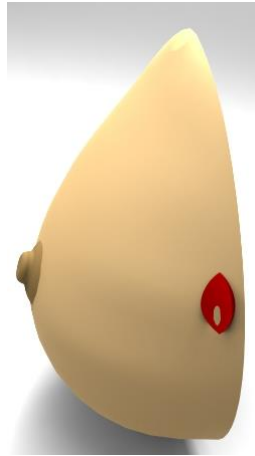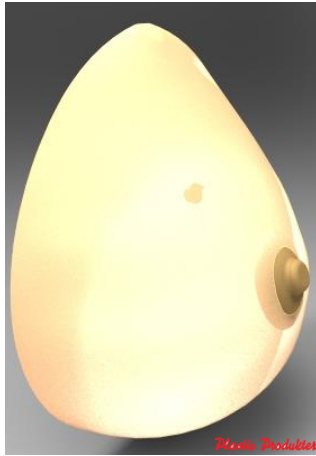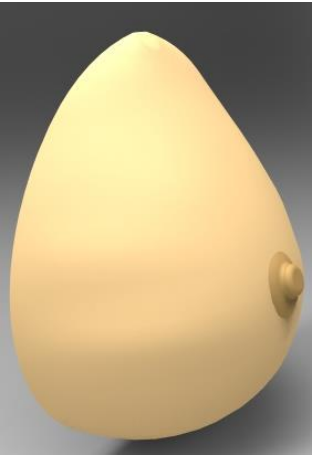

Patient has been diagnosed with unwanted tissue.  
Locate the tissue with ultrasound.  
Make a small incision in the skin.

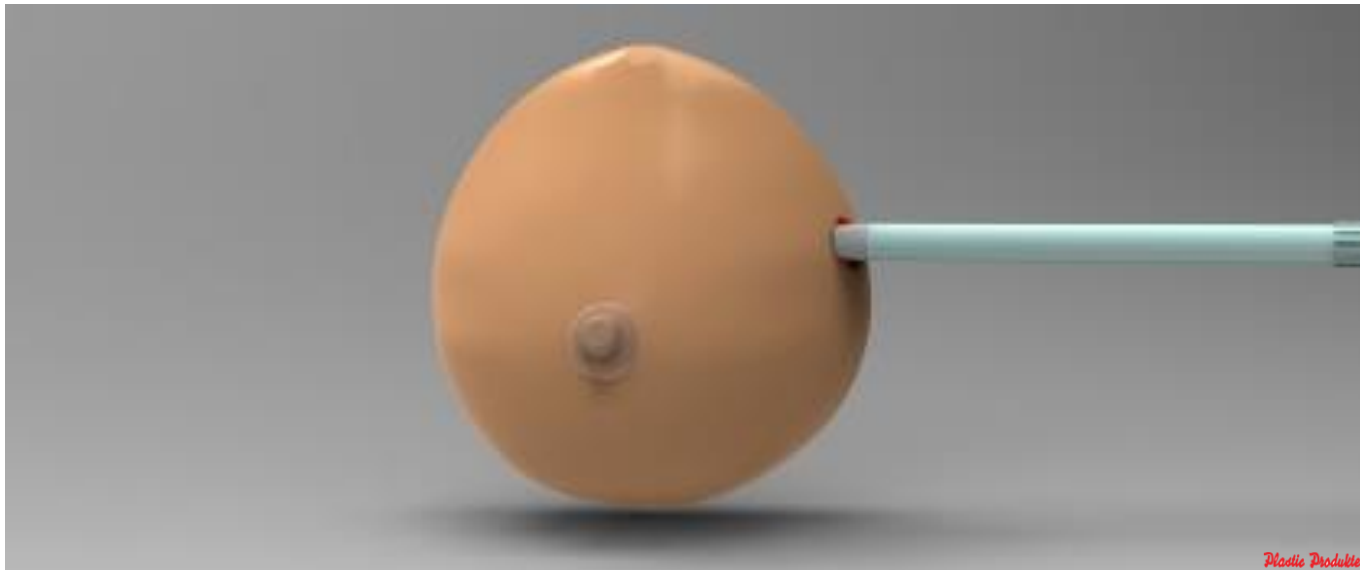


Position biopsy kit for minimally  
invasive tissue extraction.

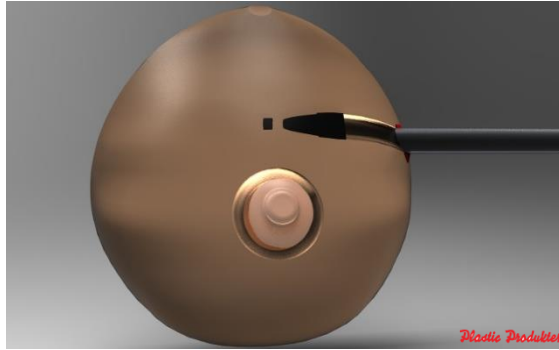


Push blunt biopsy kit towards deeply situated target tissue.



Push and rotate the sharp trephine to capture target tissue while at the same time creating vacuum.

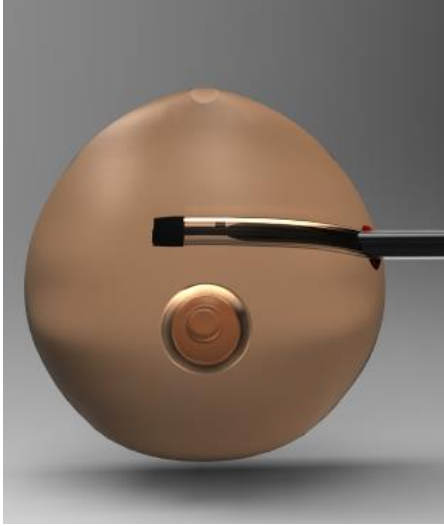

Push hooked diathermy needle through target tissue and sever it by turning the needle 360 degrees.

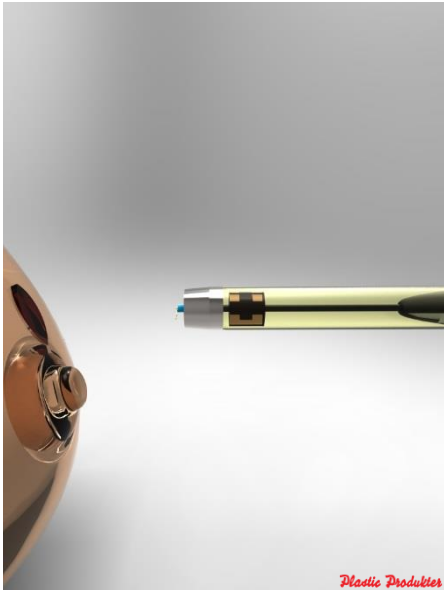

Retract biopsy kit including target tissue



# Disposable biopsy kit by NOVOAIM .

